


# ESSENTIAL WRIST BLOOD PRESSURE MONITOR



CLINICALLY  
VALIDATED

## KEY FEATURES:

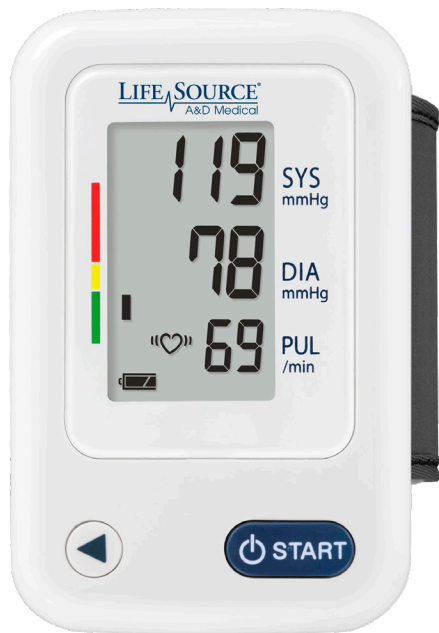

One Button Operation



AFIB Indicator



Blood Pressure Rating Indicator



60 Reading Memory



5 Year Warranty



Wrist Cuff  
13.5-21.5 cm (5.3-8.5")

A classic wrist blood pressure monitor with the benefits of speed and accuracy at a touch of a button.

## Quick and Easy Measurement

- Easy one button operation for quick and accurate measurement
- Stores 60 measurements in memory with recall
- Displays average of all measurements

## Trusted Accuracy

- Clinically validated for accuracy
- Indicates an irregular heartbeat (IHB) / Atrial fibrillation (AFib) during measurement
- Pressure Rating Indicator classifies blood pressure in accordance with World Health Organization (WHO) medically approved standards

## User Comfort

- Easy to use wrist cuff - 13.5-21.5 cm (5.3-8.5")
- Gentle inflation enables a more comfortable measurement

## Convenience

- Compact size for storage and travel

## All our blood pressure monitors come with:

- 5 year warranty
- Cuffs are not made with natural rubber latex
- Battery life indicator
- A&D Connect companion app, download for free to manually log and electronically track measurements (available on iOS and Android)

## Technical Specification

|                                     |   |
|-------------------------------------|---|
| <b>Measurement Method</b>           | Oscillometric measurement   |
| <b>Measurement Range</b>            | <b>Pressure:</b> 0-299 mmHg<br><b>Systolic pressure:</b> 60-279 mmHg<br><b>Diastolic pressure:</b> 40-200 mmHg<br><b>Pulse:</b> 40-180 beats /min |
| <b>Measurement Accuracy</b>         | <b>Pressure:</b> ±3 mmHg<br><b>Pulse:</b> ±5%   |
| <b>Power Supply</b>                 | Two (2) 1.5V alkaline batteries (LR03 or AAA)   |
| <b>Battery Life</b>                 | Approximately 250 measurements  |
| <b>Clinical Test</b>                | According to ISO81060-2:2013  |
| <b>Operating Conditions</b>         | +10 to +40 °C / 15 to 85% RH / 800 to 1060 hPa  |
| <b>Transport/Storage Conditions</b> | -20 to +60 °C / 10 to 95% RH / 700 to 1060 hPa  |
| <b>Monitor Dimensions</b>           | Approx. 56 [W] x 88 [H] x 21.5 [D] mm   |
| <b>Unit Weight</b>                  | 91 g  |
| <b>Warranty</b>                     | 5 years   |

## Packaging and Shipping Specifications

|                            | <b>Dimensions (cm) (HxLxD)</b> | <b>Weight (kg)</b> | <b>Case Cube (cubic meter)</b> | <b>UPC</b>     |
|----------------------------|--------------------------------|--------------------|--------------------------------|----------------|
| <b>Individual Unit Box</b> | 11.5 x 9.2 x 8.6               | 0.19               | 0.001                          | 093764604010   |
| <b>Case (20 units)</b>     | 25.2 x 20.2 x 45.5             | 4.05               | 0.023                          | 10093764604017 |

### What Is In The Box

- Monitor
- Instruction Manual
- Quick Start Card
- 2 AAA Batteries

**Health Canada license - 102407**

